

**Supplemental Table 1.** Contractile Properties. Percent change relative to before (5V), data are mean  $\pm$  SEM for all beats in each group CCM (10V), and After (5V).  $n = 23$ .  $*P < 0.05$ ,  $**P < 0.01$ ,  $***P < 0.001$ ,  $****P < 0.0001$ .

**Supplemental Table 2.** Calcium Handling Properties. Percent change relative to before (5V), data are mean  $\pm$  SEM for all beats in each group, CCM (10V), and After (5V).  $n = 13$ .  $*P < 0.05$ ,  $**P < 0.01$ ,  $***P < 0.001$ ,  $****P < 0.0001$ .

**Supplemental Table 3.** Electrophysiological Properties. Percent change relative to before (5V), data are mean  $\pm$  SEM for all beats in each group; CCM (10V), and After (5V).  $n = 12$ .  $*P < 0.05$ ,  $**P < 0.01$ ,  $***P < 0.001$ ,  $****P < 0.0001$ .

**Supplemental Table 4.** Baseline hiPSC-CM Properties. Absolute values of baseline (i.e., Before, 5V) cardiomyocyte excitation-contraction coupling properties for action potential (i.e., electrophysiology), calcium handling and contraction, data are mean  $\pm$  SEM;  $n = 5 - 23$ . N/A = not applicable.

**Supplemental Figure 1.** Schematic of human in vitro CCM Model. A: hiPSC-CMs are pre-plated in monolayer format on gelatin (0.1%) coated 6-well plates. B: After 2 – 28 days in culture hiPSC-CMs are dissociated and prepared for plating on Matrigel mattress substrate. C: Isolated hiPSC-CMs are plated at high density on Matrigel mattress arrayed in 48-well format (left) and assayed in [0.5 mM]  $[Ca]_0$  Tyrode solution (right). D: Commercial pulse generator and standard clinical CCM pulse parameters (right) are used stimulate hiPSC-CMs, cardiac function is assessed by video and florescence measurements (left). E: Representative contraction recording before CCM (5V), CCM (10V) and after (5V).

**Supplemental Figure 2.** Percent Change for the Effect of Extracellular Calcium Modulation on CCM Response. hiPSC-CMs were exposed to increasing concentrations of extracellular Ca [ $\text{Ca}_0$ ] 0.25 – 2 mM. Summary bar graphs of immediate effects.  $n = 6 - 8$  per group. Transformed data from figure 5.  $*P < 0.05$ ,  $**P < 0.01$ ,  $***P < 0.001$ ,  $****P < 0.0001$ .

**Supplemental Figure 3.** Electric Field Numerical Modeling. A: Prospective and top view of the geometry of the platinum electrodes inserted in a glass bottom well. B: Electric field intensity in the YZ plane perpendicular to the electrodes for 1 V applied. For quantitative analysis, the values of the  $|E|$  field were extracted in the region of interest (ROI). C:  $|E|$  field along the y axis; the ROI is highlighted by a rectangle. D: Scaling table to convert the  $|E|$  field in the ROI obtained at 1 V applied to significant experimental values, 5 V and 10 V, respectively.

**Supplemental Figure 4.** Effect of Voltage Amplitude on CCM Response. A: Contraction traces for each group, Before (5V), CCM (0 - 15V), After (5V), hiPSC-CMs were exposed to increasing CCM pulse amplitude (0 – 15) Volts. [ $\text{Ca}_0$ ] 0.5 mM. B: Transformed data demonstrating the effect of CCM pulse amplitude on contraction amplitude During CCM (Hill slope = 8.0).  $n = 1$ .  $\text{EV}_{50}$  = Effective Voltage 50%.

**Supplemental Figure 5.** Effect of 0V CCM on hiPSC-CM Contractility. A: Contraction traces for each group, Before (5V), CCM (0V), After (5V). [ $\text{Ca}_0$ ] 0.5 mM.

**Supplemental Video.** hiPSC-CMs on Matrigel Mattress. Matrigel Mattress in one well of 48-well glass bottom plate 4x. hiPSC-CMs form monolayer morphology and robust contraction at ~ 24 hours post plating. White arrow indicated edge of Matrigel mattress. Scale bar, 1 mm.

**Supplemental Table 1. Contractile Properties**

<b>Parameter</b>	<b>CCM</b>	<b>After</b>
<b>Amplitude</b>	<b>16 ± 4%**</b>	<b>4 ± 5%</b>
<b>Time to Peak 50%</b>	<b>-20 ± 9%*</b>	<b>7 ± 5%</b>
<b>Time to Peak 90%</b>	<b>-22 ± 8%*</b>	<b>6 ± 5%</b>
<b>Time to Baseline 50%</b>	<b>-8 ± 5%</b>	<b>4 ± 4%</b>
<b>Time to Baseline 90%</b>	<b>-12 ± 6%*</b>	<b>5 ± 5%</b>
<b>Contraction Duration 10%</b>	<b>-13 ± 6%</b>	<b>3 ± 5%</b>
<b>Contraction Duration 50%</b>	<b>-6 ± 5 %</b>	<b>3 ± 5%</b>
<b>Contraction Duration 90%</b>	<b>0 ± 5%</b>	<b>3 ± 4%</b>
<b>N</b>	<b>23</b>	<b>23</b>

**Supplemental Table 2. Calcium Handling Properties**

<b>Parameter</b>	<b>CCM</b>	<b>After</b>
<b>Amplitude</b>	<b><math>13 \pm 5\%^*</math></b>	<b><math>-10 \pm 2\%^{**}</math></b>
<b>Time to peak</b>	<b><math>-22 \pm 3\%^{****}</math></b>	<b><math>-1 \pm 3\%</math></b>
<b>Ca Rise Time</b>	<b><math>-33 \pm 3\%^{****}</math></b>	<b><math>5 \pm 2\%^*</math></b>
<b>Ca Duration 50%</b>	<b><math>-10 \pm 2\%^{***}</math></b>	<b><math>0 \pm 1\%</math></b>
<b>Ca Duration 90%</b>	<b><math>-2 \pm 1\%</math></b>	<b><math>1 \pm 1\%</math></b>
<b>N</b>	<b>13</b>	<b>13</b>


**Supplemental Table 3. Electrophysiological Properties**

<b>Parameter</b>	<b>CCM</b>	<b>After</b>
<b>TRise</b>	<b>-13 ± 5%*</b>	<b>-6 ± 17%</b>
<b>APD50</b>	<b>-8 ± 2%**</b>	<b>18 ± 6%</b>
<b>APD75</b>	<b>-9 ± 1%****</b>	<b>18 ± 4%</b>
<b>APD90</b>	<b>-10 ± 1%****</b>	<b>17 ± 5%</b>
<b>n</b>	<b>12</b>	<b>12</b>

**Supplemental Table 4. Baseline hiPSC-CM Properties**


Parameter	Action Potential	Calcium	Contraction
Amplitude (a.u.)	N/A	0.20 ± 0.03	1882.5 ± 187.2
Time to Peak 50% (ms)	N/A	N/A	166.4 ± 15.5
Time to Peak 90% (ms)	N/A	272.6 ± 15.0 <sup>#</sup>	245.4 ± 18.1
TRise (ms)	70.0 ± 23.9	100.5 ± 7.2	N/A
Contraction / CaT (Up), Slope (a.u./s)	N/A	22104 ± 5517	20088.6 ± 1920.0
Time to Baseline 50% (ms)	N/A	N/A	666.2 ± 29.0
Time to Baseline 90% (ms)	N/A	N/A	807.5 ± 21.1
Relaxation or CaT (Down), Slope (a.u./s)	NA	-8926 ± 2144	-9055.8 ± 843.8
Duration 10% (ms)	N/A	N/A	375.4 ± 21.1
Duration 50% (ms)	519.7 ± 24.2	517.1 ± 17.0	500.0 ± 19.3
Duration 75% (ms)	581.1 ± 24.2	N/A	N/A
Duration 90% (ms)	613.7 ± 19.4	676.0 ± 10.3	561.8 ± 13.0
Beat Rate, Spontaneous (BPM)	N/A	N/A	33.3 ± 6.4
Interval, Spontaneous (ms)	N/A	N/A	2353.2 ± 472.4
Interval, Paced 1 Hz (ms)	N/A	N/A	998.18 ± 2.7
Beat Rate, Paced 1 Hz (BPM)	60	60	60

<sup>#</sup> = Time to 100%

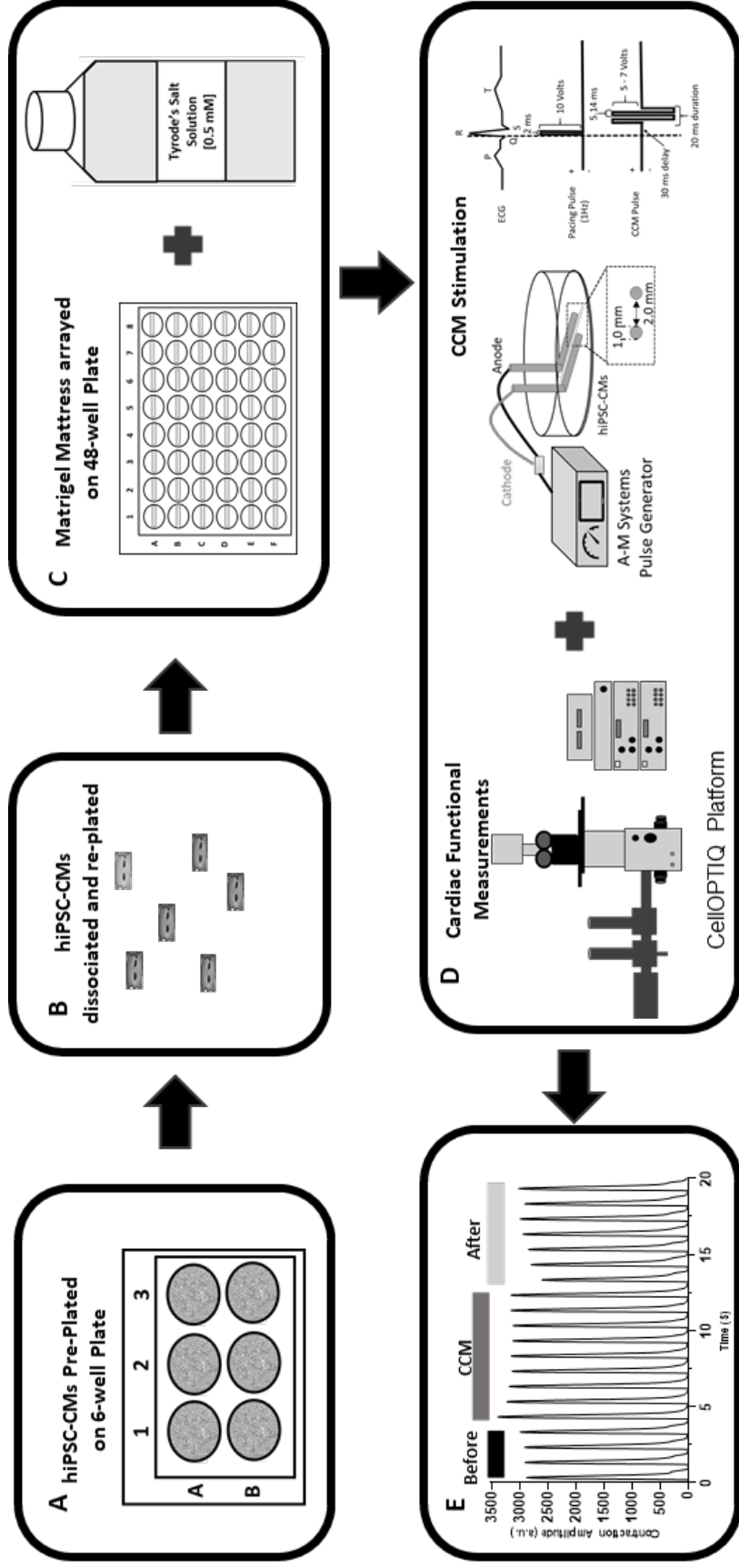


Supplemental Figure 1


2  $\text{Ca}_0$  [mM]

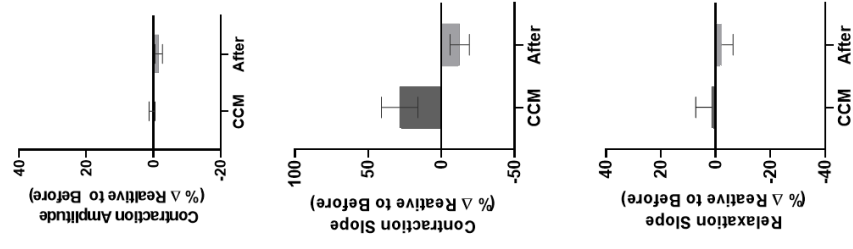


1

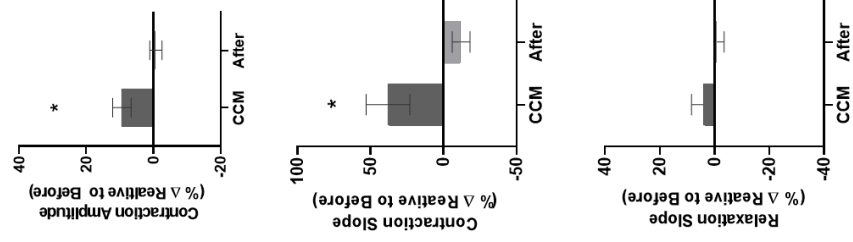

0.5

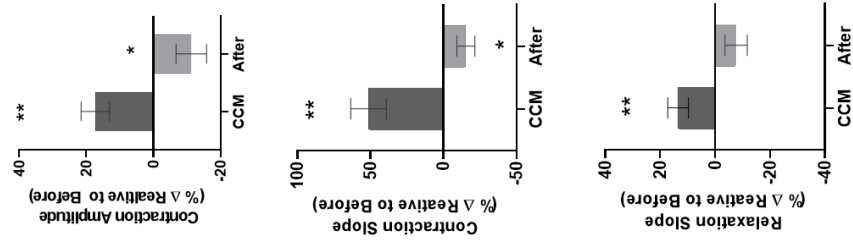


0.25

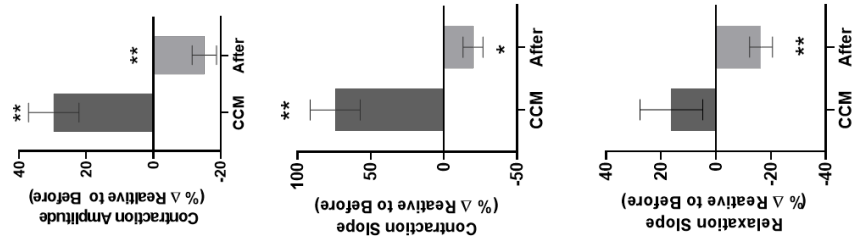

Supplemental Figure 2




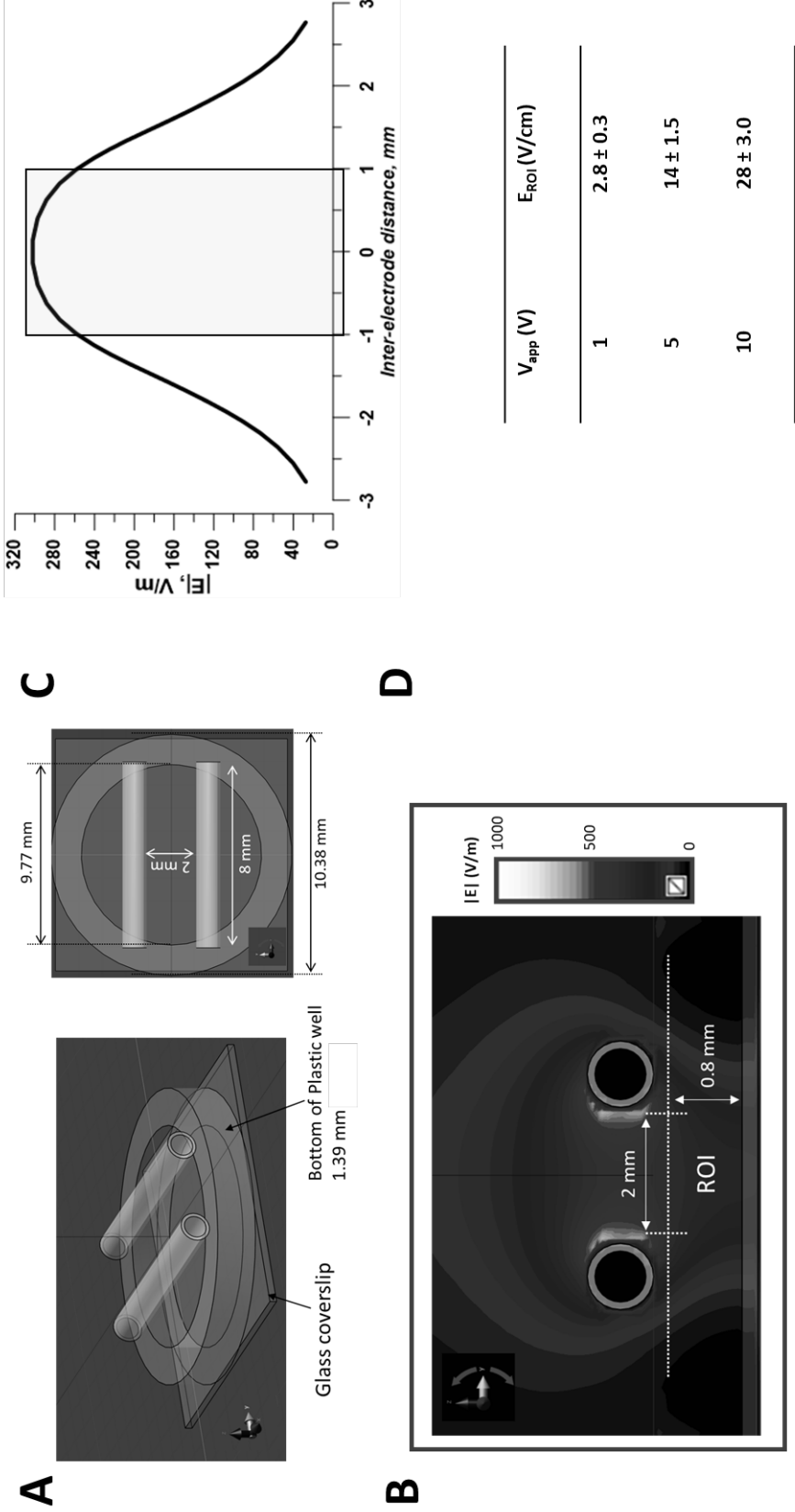


Supplemental Figure 3

A


CCM Amplitude Voltage

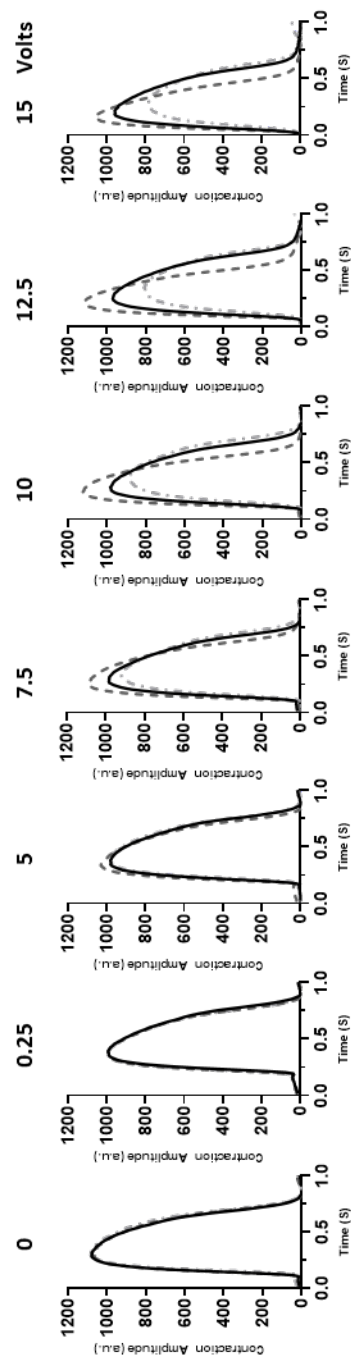


B

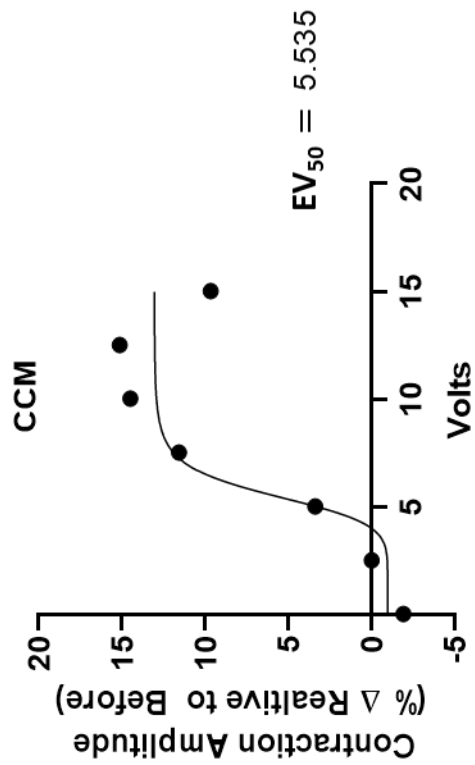

Supplemental Figure 4

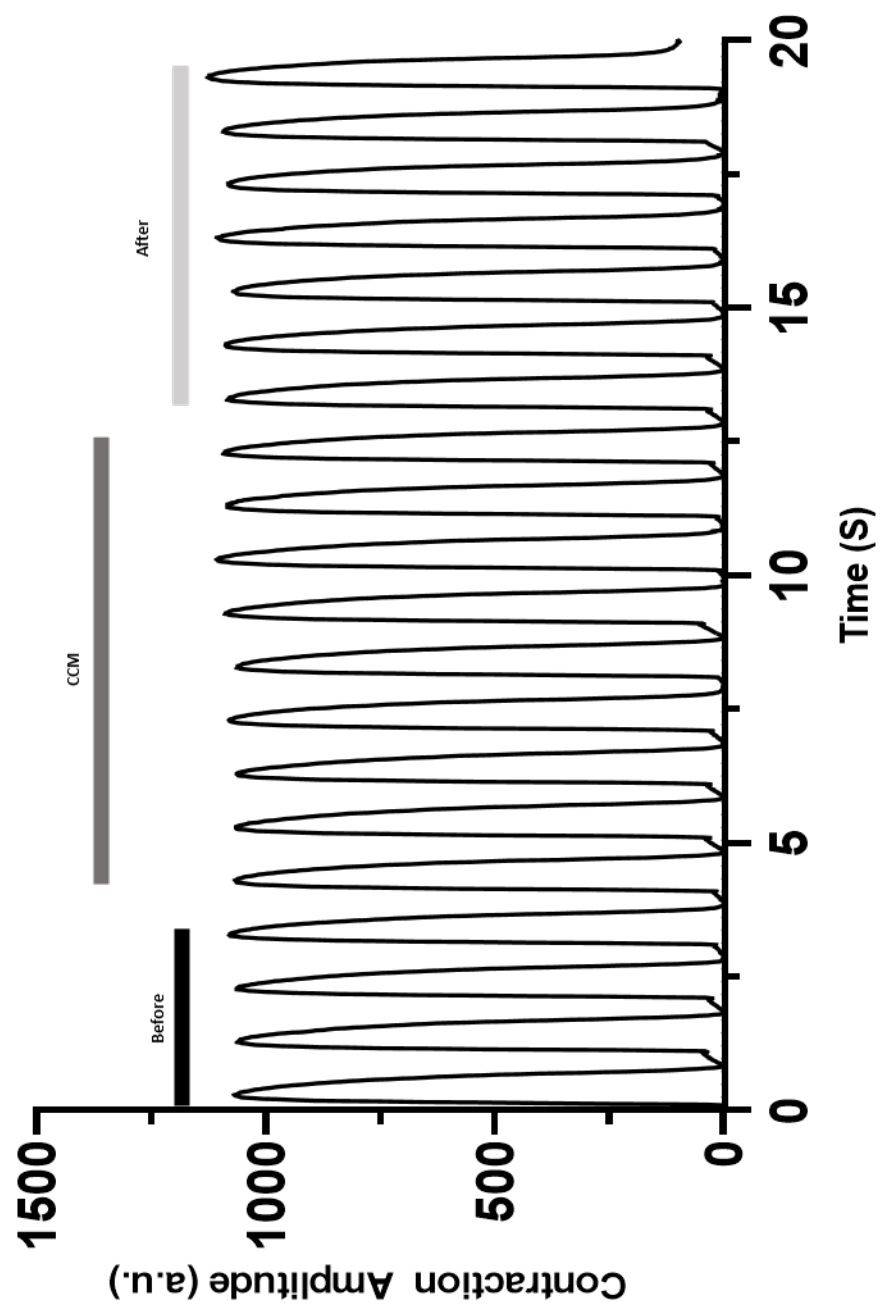

Supplemental Figure 5

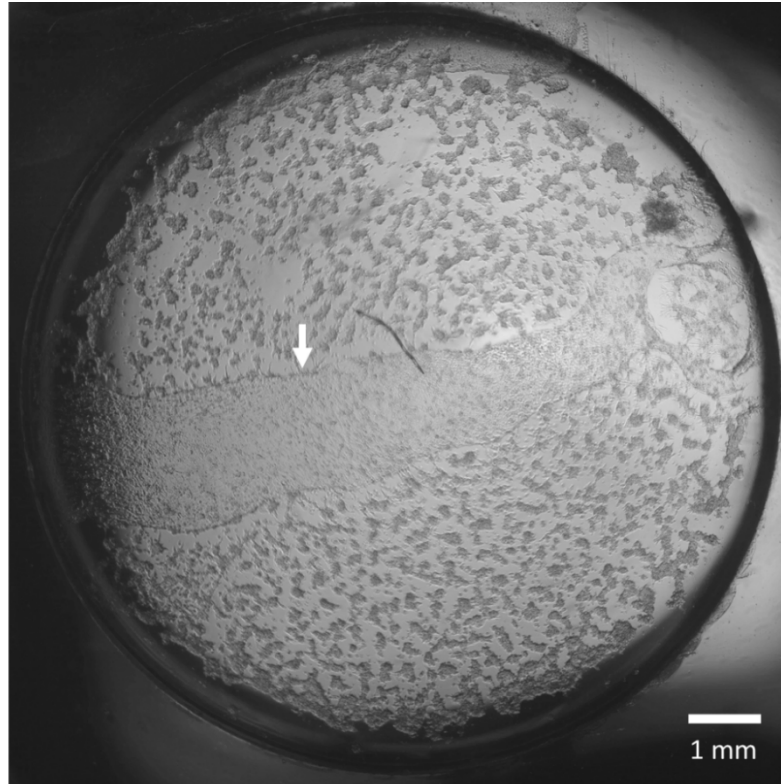

**Supplemental Video**